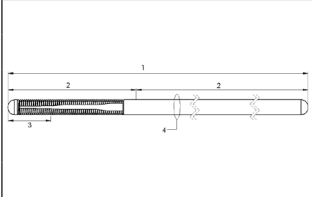
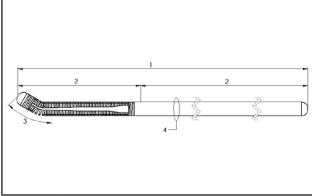
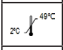
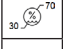
# Champion

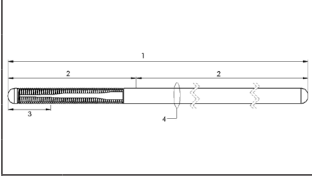
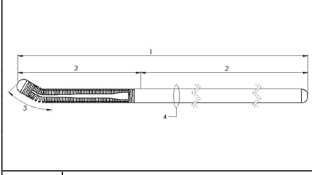
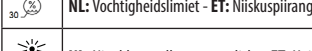
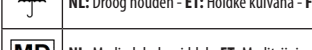
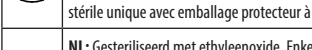
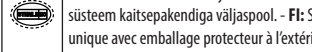
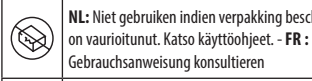
## PTCA Guide Wire

EN	<i>English</i>	<b>PTCA Guide Wire</b>
BG	<i>Bulgarian/Български език</i>	<b>Водач за PTCA</b>
HR	<i>Croatian/Hrvatski</i>	<b>Žica vodilica PTCA</b>
CS	<i>Czech/Česky</i>	<b>Vodící drát PTCA</b>
DA	<i>Danish/Dansk</i>	<b>PTCA guidewire</b>
NL	<i>Dutch/Nederlands</i>	<b>PTCA voerdraad</b>
ET	<i>Estonian/Eesti</i>	<b>PTCA juhtetraat</b>
FI	<i>Finnish/Suomi</i>	<b>PTCA ohjainlanka</b>
FR	<i>French/Français</i>	<b>Fil-guide ACTP</b>
DE	<i>German/Deutsch</i>	<b>PTCA Führungsdraht</b>
EL	<i>Greek/Ελληνικά</i>	<b>Οδηγό σύρμα PTCA</b>
HU	<i>Hungarian/Magyar</i>	<b>PTCA vezetődrót</b>
IT	<i>Italian/Italiano</i>	<b>Filo guida PTCA</b>
LV	<i>Latvian/Latviešu</i>	<b>PTCA vadītājstīga</b>
LT	<i>Lithuanian/Lietuvių k.</i>	<b>PTCA kreipiamoji viela</b>
NO	<i>Norwegian/Norsk</i>	<b>PTCA ledesonde</b>
PL	<i>Polish/Polski</i>	<b>Prowadnik PTCA</b>
PT	<i>Portuguese/Potuguês</i>	<b>Fio-guia de PTCA</b>
RO	<i>Romanian/Română</i>	<b>Fir de ghidaj PTCA</b>
SK	<i>Slovak/Slovenčina</i>	<b>Vodiaci drôt PTCA</b>
SL	<i>Slovenian/Slovenščina</i>	<b>Vodilna žica PTCA</b>
ES	<i>Spanish/Español</i>	<b>Alambre guía para ACTP</b>
SV	<i>Swedish/Svenska</i>	<b>PTCA Ledare</b>
TR	<i>Turkish/Türkçe</i>	<b>PTCA Kılavuz Tel</b>

**EN: Explanation of symbols used on the package labels - BG: Значение на символите върху етикетите на опаковката - HR: Objašnjenja simbola koji se koriste na naljepnici na pakiranju - CS: Vysvětlení symbolů použitých na štítcích obalů - DA: Forklaring af symboler anvendt på emballagens etiketter**

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## Champion PTCA Guide Wire

### Description:

The Champion PTCA Guide Wire is characterized by a stainless-steel shaft with PTFE coating and a 3 cm platinum tip with hydrophilic coating. The guide wire is EO sterilized, non-pyrogenic and available in the following variants:

Subtype	Tip subtype	Tip stiffness	Stent zone stiffness	Available variants
CAR	<ul style="list-style-type: none"> <li>• Straight</li> <li>• J-tip</li> </ul>	<ul style="list-style-type: none"> <li>• Extra floppy</li> <li>• Floppy</li> <li>• Intermediate</li> <li>• Standard</li> </ul>	<ul style="list-style-type: none"> <li>• Standard</li> <li>• Stiff</li> </ul>	<ul style="list-style-type: none"> <li>• Tip shape</li> <li>• Tip stiffness</li> <li>• Guide wire length</li> </ul>
CTO	<ul style="list-style-type: none"> <li>• Straight</li> <li>• J-tip</li> </ul>	<ul style="list-style-type: none"> <li>• 3 gf</li> <li>• 6 gf</li> <li>• 12 gf</li> </ul>	<ul style="list-style-type: none"> <li>• Standard</li> </ul>	

For further information, see SP Medical's website.

### Intended purpose:

The Champion PTCA Guide Wire is indicated for general use in the coronary vasculature to aid in the selective placement of interventional devices during diagnostic and/or therapeutic procedures.

### Contraindications:

**Not for use in the cerebral vasculature.**

### Warnings:

1. SINGLE USE. This guide wire is intended for single use only. Reusing the guide wire involves high risk of contamination and locking of the guide wire inside the interventional device or human body due to guide wire kinking or wear-out of the hydrophilic coating.
2. This guide wire is provided in sterile condition. Do not re-sterilize or re-use the guide wire.
3. Do not withdraw through a metal needle.
4. Avoid interaction with metal / hard plastic devices designed with sharp edges or burrs that may separate the coating from the guide wire.
5. Never twist or apply excessive force as the core may penetrate the coil and cause vessel damage.
6. Never advance, withdraw or torque the guide wire against resistance until the cause of the resistance is determined by fluoroscopy.
7. Do not attempt to move the guide wire without observing the resultant tip response.
8. Do not use alcohol, antiseptic solutions, other solvents or dry gauze on the guide wire, as this may damage the hydrophilic coating.
9. The guide wire contains a metallic core. Do not use with any inappropriate equipment (e.g., MRI).
10. To be used before the expiry date stated on the package.

### Precautions:

1. Prior to opening, the sterile package should be carefully inspected to see if it is still intact. The inspection must be carried out just before the guide wire is used and in good light conditions. Sealing quality on the transparent side of the peel pouch must be carefully inspected for channels and cracks. Both front and back side must be inspected for pinholes or other packaging defects such as abrasions or cracks. Use a magnification tool if in doubt of sterile barrier integrity.
2. Do not use if the package / sterile barrier is broken / damaged.
3. Prior to use, carefully inspect the guide wire for damages.
4. Do not use damaged guide wires.
5. Champion PTCA Guide Wire for highly stenosed lesions (CTO) is designed with a stiff tip. Use lowest tip stiffness that will treat the lesion and take due care to minimize the risk of perforation or damage to blood vessels.
6. SP Medical does not recommend a particular technique or device configuration for the use, and the steps contained in these instructions for use are for information purposes only.
7. The guide wire should only be selected and used by experienced physicians, trained in invasive techniques, the use of guide wires and familiar with side effects and hazards commonly associated with interventional procedures.
8. The physician should evaluate the appropriateness according to the individual patient condition and his or her medical training and experience.
9. It is always the physician's responsibility to determine and ensure the patient's suitability for the percutaneous coronary intervention (PCI) where the guide wire is used.

### Potential side effects:

Possible complications include but are not limited to the following:

- Vessel / tract trauma
- Vessel / tract narrowing / occlusion
- Infection and inflammation
- Allergic reaction

### Serious incident:

Every serious incident that has occurred in relation to the device should be reported without delay to the manufacturer and the competent authority of the Member State in which the user and / or patient is established. For further information, please contact the manufacturer.

### Compatibility:

The guide wire may be used in combination with interventional devices in a CathLab environment. The end-user is solely responsible for selection of the appropriate device. There should be at least 0.0004" (0.01 mm) clearance between the lumen of the interventional device and the guide wire.

### Preparations before use:

1. Carefully remove the guide wire from the dispenser by releasing the wire from the retention clip and pushing the guide wire in the open area of the dispenser.
2. Inspect the guide wire thoroughly to make certain that it is not kinked or otherwise damaged.
3. Confirm the compatibility of the guide wire diameter with the interventional device before use.

### Directions for use:

1. If desired, shape the guide wire tip using standard technique. When preparing to shape the tip section, hold the tip at the middle joint and between two fingers. Gently "brush" the spring coil to identify the plane of flexure and complete the shaping procedure in this plane.
2. Moisten the guide wire with sterile saline to increase the surface lubricity.
3. Introduce the guide wire using a guide wire introducer.
4. Remember to remove the guide wire introducer by sliding it over the proximal end of the guide wire. Secure the guide wire within the haemostatic valve, being careful not to over-tighten the compression fitting.
5. A torque application device may be applied to the proximal end of the guide wire, if needed.
6. Between uses during the same procedure, place the guide wire in a saline-filled container to keep it moist. Avoid wiping with damp cloths as particles from the cloth can adhere to the surface of the guide wire.
7. When the procedure is finished, remove the guide wire and interventional device according to procedural protocol.
8. Discard the product after use according to local instructions for hazardous waste.

### If using a rapid-exchange system:

1. When using a rapid-exchange balloon catheter on the guide wire, advance the guide wire into the position according to manufacturers' instructions.
2. Manifold flushing of the catheter lumen during guide wire advancement is recommended to remove any residual air.

### Storage:

To be stored under cool, dark and dry conditions.  
Humidity limits: 30% - 70%. Temperature limits: 2°C - 49°C

### Expiry date:

The guide wire must be used no later than stated on the label.

### Disposal:

Disposal shall be handled in accordance with accepted medical practice and applicable local laws and regulations for hazardous waste.

### Safety information:

The IFU is delivered in hard copy together with the products and kept updated and available for electronic download on SP Medical's website.

Summary of safety and clinical performance can be obtained via Eudamed using SP Medical as legal manufacturer and the basic UDI-DI: 5710086-GWC2J2.

<https://ec.europa.eu/tools/eudamed/#/screen/search-device>

### Disclaimer of liability:

SP Medical A/S is not liable for defects / deterioration resulting from abnormal use or modifications made to the product. Under these circumstances, the product is not covered by the guarantee. SP Medical A/S disclaims liability for direct or indirect injuries that may occur as a consequence of the product being modified or wrongly used.

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